

# MV Set for Blood Transfusion (Hemoset)



Bulb Latex



Automatic shut off  
floating valve



Smooth Flow  
Controller

## MV Set for Blood Transfusion

The Measured Volume Set (M.V Set) is Cylindrical in shape, which has a flexible, clear and calibrated measured volume chamber. Its chamber has an in built Injection port and air inlet with bacteria filter.

---

Automatic shut off floating valve


---

Smooth Flow Controller

---

Transparent, cylindrical & calibrated chamber

---

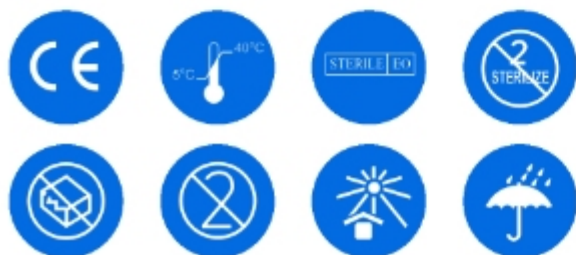

## Features

It have an automatic shut off floating valve provided in chamber to prevent air trapping and to indicate volume of blood in the chamber.

- Measured volume chamber in size 100/150 ml
- 15 micron fluid filter in the drip chamber
- 20 Drops equivalent to 1ml  $\pm$  0.1 ml (1g  $\pm$  0.1g)