
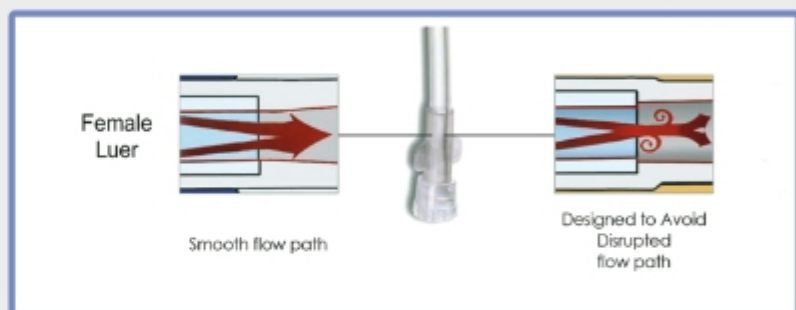
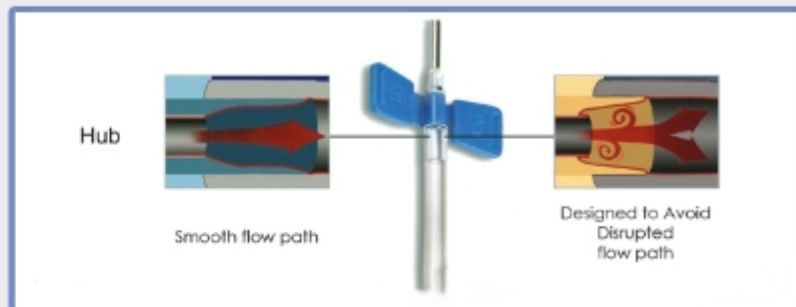


AV FISTULA NEEDLE

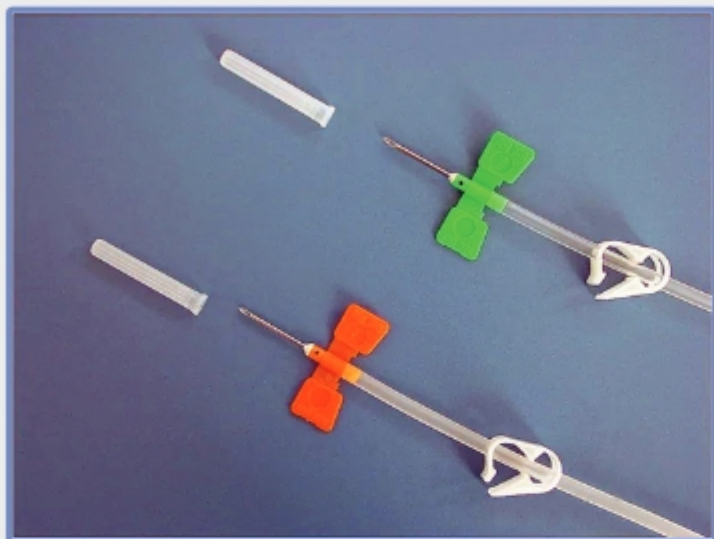

- When flipping the needle is necessary, orientation dots on the hub allow you to aim the bevel precisely and reduce the risk of tissue damage.
- No two patients are the same: GMS offers you choice of gauges, and tubing lengths.

Features

- Ultra-thin wall stainless steel needle (with or without back-eye) siliconized on both internal and external surfaces.
- Provides optimum penetration profile.
- Minimizes pain during insertion , Maximizes smooth, laminar flow
- Soft, kink-resistant tubing made from medical-grade PVC.
- Color-coded wings facilitate flipping of needle while allowing for easy gripping and fixation.
- Orientation dot on needle hub enables precise maneuvering of bevel.
- Luer has universal 6% taper to allow smooth connection with standard bloodline sets.
- Available in single pack or twin-pack for better use of space.
- Available with Red, Blue or White Clamp.
- Available from 15 cm to 30 cm Tubing.
- Sizes : 15 G, 16 G & 17 G.